



COOLCELL





NEW

HIGH RISK




Features & Benefits

- Cooling additive inside Gel, helps draw heat away from the body
- Unique combination of AIR and Cooling GEL technology effectively provides pressure relief, comfort and cooling to the high risk user
- 34 interconnected air-filled cells allow movement of air across surface of cushion, each covered with a layer of pressure relieving silicone-based Cooling Gel
- Air cushion enhances capillary function
- **Silicone-based TruGel**
 - enhances pressure re-distribution
 - minimizes the effects of shear forces
 - has a positive effect on skin temperature
- When pressure is applied, air circulates to each of the cells, distributing weight evenly, thus reducing high pressure concentrations around bony prominences
- The Coolcell is made up of 3 layers of polyurethane (pu) film giving a total thickness of 1mm approx making it highly durable and difficult to puncture
- Clear protective pu film coating is easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Easy to inflate and adjust
- Includes outer cover and pump

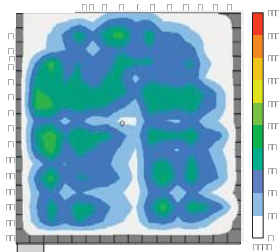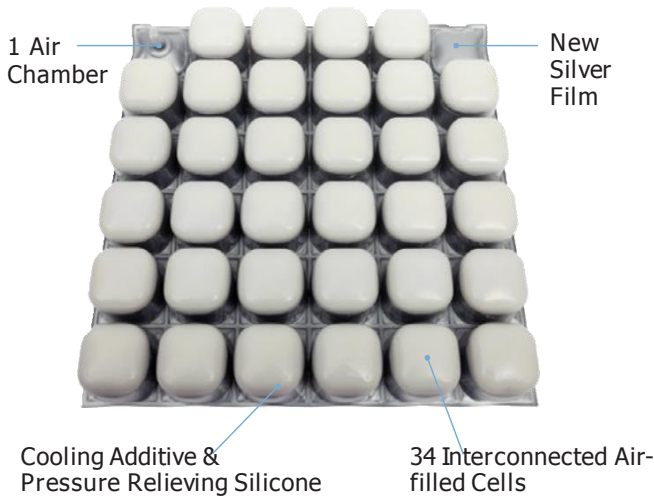

Pressure diagram with RELAX CoolCell cushion

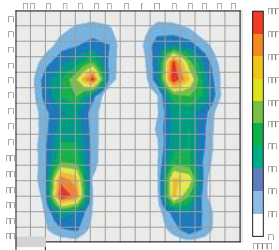

Pressure diagram without wheelchair cushion


CoolCell cushion with spacer cover

Product Code	Dimensions (Width x Length x Height)	Cushion Weight
RC002-1616	40cm x 40cm x 5cm/ 16" x 16" x 2"	1.5kg/3.3lb
RC002-1717	43cm x 43cm x 5cm/ 17" x 17" x 2"	1.6kg/3.5lb
RC002-1816	45cm x 40cm x 5cm/ 18" x 16" x 2"	1.6kg/3.5lb
RC002-1818	45cm x 45cm x 5cm/ 18" x 18" x 2"	1.7kg/3.7lb

CoolCell cushion with standard black outer cover

Product Code	Dimensions (Width x Length x Height)	Cushion Weight
RC004-1616	40cm x 40cm x 5cm/ 16" x 16" x 2"	1.5kg/3.3lb
RC004-1717	43cm x 43cm x 5cm/ 17" x 17" x 2"	1.6kg/3.5lb
RC004-1816	45cm x 40cm x 5cm/ 18" x 16" x 2"	1.6kg/3.5lb
RC004-1818	45cm x 45cm x 5cm/ 18" x 18" x 2"	1.7kg/3.7lb

The weight limit for the CoolCell cushion is 130kgs/287lbs

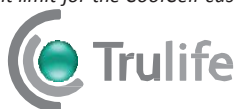



Our Relax range of wheelchair cushion covers are all designed to offer protection to the user and the product. The covers are designed to be fluid repellent, vapour permeable, stretchable and hand-washable.

CoolCell

Spacer Fabric Cover



Materials:

Topside: Spacer fabric material enhances redistribution of pressure peaks and dispersion of air and water.

Underside: Weft knitted nylon polyurethane griplock, treated against flammability, protects against fungi and bacteria.

Features:

Excellent dispersion of air and moisture
Enhances redistribution of pressure peaks
Millions of microfilaments act as miniature pressure springs
Highly durable
Washable at 60°C
Velcro included to prevent migration in wheelchair
High grip/anti-slip surface

Codes & Sizes

CC002-1616
CC002-1717
CC002-1816
CC002-1818



Ability in Motion



1800 994 408
Fax 5335 1604
admin@abilityinmotion.com.au