





HIGH RISK

Features & Benefits


- Two separate **air chambers** provide greater stability and pressure redistribution for buttocks and Ischial Tuberosities
- Air cushion enhances capillary function
- **Silicone-based TruGel**
 - enhances pressure re-distribution
 - minimizes the effects of shear forces
 - has a positive effect on skin temperature
- The cells are **contoured** towards the rear making the G2 the perfect all round pressure relieving and positioning cushion
- When pressure is applied, air circulates to each cell in both chambers distributing weight evenly. This process reduces high pressure concentrations around bony prominences, preventing skin breakdown
- The Gelcell G2 is made up of **3 layers of polyurethane** (pu) film giving a total thickness of 1mm approx making it highly durable and difficult to puncture
- Clear protective pu film coating is easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Easy to inflate and adjust
- Includes outer cover and pump

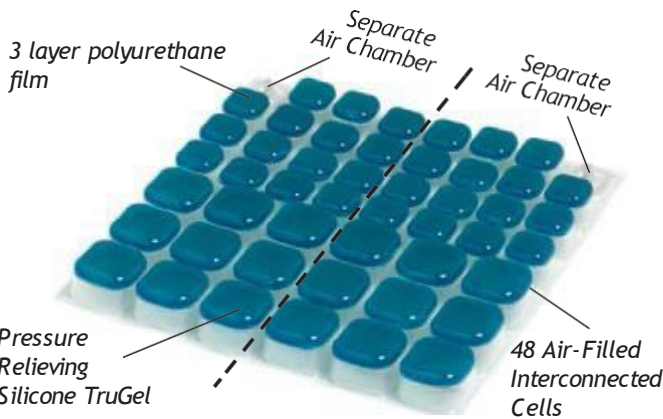

Contoured Design at Rear of Cushion

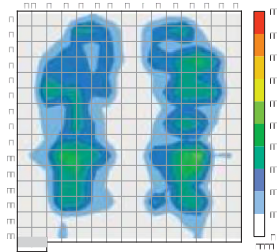

Pressure diagram with RELAX Gelcell G2 cushion

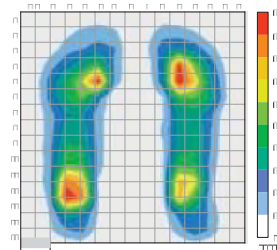

Pressure diagram without wheelchair cushion

Product Code	Dimensions (Width x Length x Height)	Cushion Weight
RG005-1414	36cm x 36cm x 4cm/14"x14"x1.5"	1.3kg/2.9lb
RG005-1515	38cm x 38cm x 4cm/15"x15"x1.5"	1.3kg/2.9lb
RG005-1616	40cm x 40cm x 4cm/16"x16"x1.5"	1.5kg/3.3lb
RG005-1618	40cm x 45cm x 4cm/16"x18"x1.5"	1.6kg/3.5lb
RG005-1717	43cm x 43cm x 4cm/17"x17"x1.5"	1.6kg/3.5lb
RG005-1816	45cm x 40cm x 4cm/18"x16"x1.5"	1.6kg/3.5lb
RG005-1818	45cm x 45cm x 4cm/18"x18"x1.5"	1.6kg/3.5lb

The weight limit for the Gelcell G2cushion is 130kgs/287lbs



1800 994 408

Fax 5335 1604

admin@abilityinmotion.com.au